Solanas[®] Posterior Cervico-Thoracic Fixation System

۲

SYSTEM FEATURES

- Top loading, top tightening system
- Multi-axial screws offer up to 80 degree variability
- Proven buttress-thread closure mechanism helps prevent cross-threading

۲

- User friendly, color-coded implants and instrumentation ease identification
- Integrated functioning of screws, rods, and hooks





CERVICAL | Data Sheet

۲

The precision-machined size and pitch of the locking cap threads are designed to assure quick insertion while minimizing potential for cross threading.



Adjustable Bridges in lengths 26, 37 & 60mm connect rods to increase construct stability.


۲

Rod-to-Rod Connectors offer additional system variability by adjoining the 3.3mm rod with the 5.5mm rod of the Zodiac[®] system.



Implants (in millimeters)

۲

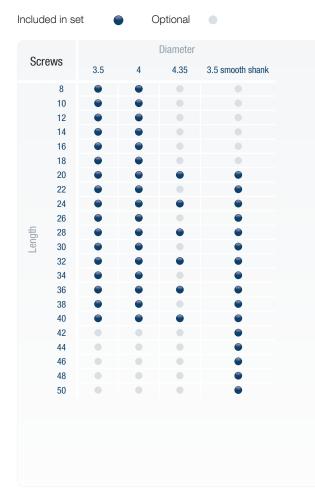

Part No.	Description	Qty.
Rods		
63001-120 63001-240 63002-500	Solanas Rod - 120mm Solanas Rod - 240mm 3.3mm - 5.5mm Transitional Rod - 500mm	4 4 4
Adjustable Bridges		
63003-26 63003-37 63003-60	26mm Adjustable Bridge - Small 37mm Adjustable Bridge - medium 60mm Adjustable Bridge - Large	2 2 2
Set Screws		
63020	Solanas Hexalobe Set Screw	20
Rod-to-Ro 63008 63009 63010 63011	od Connectors 3.3mm - 5.5mm Axial Rod Connector 3.3mm - 5.5mm Parallel Rod Connector 3.3mm - 3.3mm Parallel Rod Connector 3.3mm - 3.3mm Axial Rod Connector	2 2 2 2
Cervical Hooks		
63012 63013	Cervical Hook - 6mm Cervical Hook - 5mm	6 6
Lateral Connectors		
63014 63015	Closed Lateral Rod Connector Open Lateral Rod Connector	6 6

X Alphatec Spine°

CORPORATE HEADQUARTERS 5818 El Camino Real Carlsbad, California 92008

۲

CUSTOMER SERVICE

Toll Free: 800.922.1356 Local: 760.431.9286 Fax: 800.431.1624 Solanas, Zodiac and the Alphatec Spine logo are trademarks of Alphatec Spine, Inc. © 2015 Alphatec Spine, Inc. All rights reserved.

LIT-63101G | 2015/07/10

ReBrand_LIT-63101G_Solanas_DS_2p_071015.indd 2

7/10/2015 3:20:21 PM

۲